



STARMed Ventumask 30

For CPAP therapy with Venturi flow driver, adjustable PEEP valve and pressure manometer



Now available

with nebuliser port option

CV205S (Small)

CV205M/L (Medium/Large)

CV205XL (Extra Large)

Nebuliser kit not included

Critical Care ▪ Patient Interfaces ▪ STARMed


STARMed Ventumask 30 with Venturi flow driver & adjustable PEEP valve

StarMed Ventumask 30

StarMed Ventumask 30 is a unique design for CPAP therapy that is simply activated by connecting to an appropriate oxygen supply. It can be used in a hospital environment or pre-hospital use, for example in an ambulance.

Advantages and Technical Features

- Ventumask 30 is a complete and easy to use system, where the oxygen supply tubes, a Venturi flow driver and a PEEP valve are integral to a CPAP mask fitted with a pressure gauge
- The built-in adjustable PEEP valve can provide patient pressures in a range of settings up to 20 cmH₂O
- The integrated Venturi flow driver can generate patient flows up to 80 L/min while consuming a reduced amount of compressed oxygen
- O₂% can be adjusted from 30% to 100% at each PEEP level
- The user is guided by a performance table integral to the product, for independent PEEP and O₂% settings

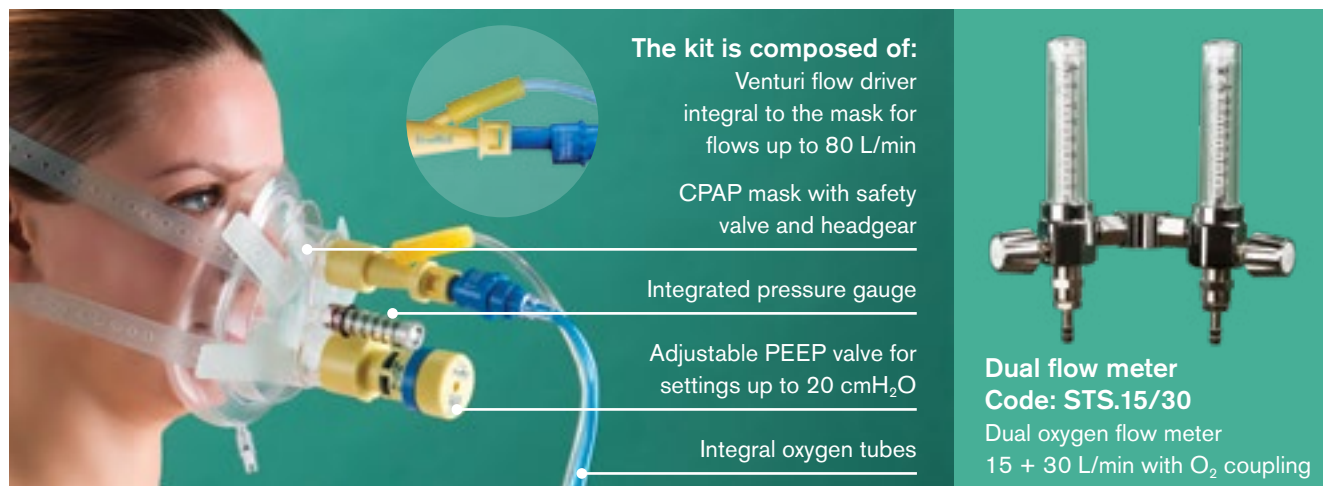

Ventumask 30 for CPAP therapy with Venturi flow driver, adjustable PEEP valve and pressure manometer

Requirement		Setting	Flow delivered
O ₂ % approx	PEEP cmH ₂ O	Flow O ₂ L/min A + B	Flow approx L/min
30%	5	12 + 0	78
	7.5	12 + 0	75
	10	12 + 0	70
	12.5	12 + 0	65
	15	12 + 0	60
40%	20	14 + 0	60
	5	11 + 8	78
	7.5	11 + 8	72
	10	12 + 6	70
	12.5	12 + 4	65
	15	12 + 4	60
	20	13 + 2	60

Requirement		Setting	Flow delivered
O ₂ % approx	PEEP cmH ₂ O	Flow O ₂ L/min A + B	Flow approx L/min
50%	5	10 + 18	80
	7.5	10 + 18	77
	10	10 + 16	72
	12.5	10 + 14	65
	15	10 + 12	60
60%	20	12 + 10	60
	5	9 + 26	80
	7.5	9 + 26	74
	10	10 + 24	75
	12.5	10 + 22	72
	15	10 + 22	70
	20	10 + 18	60
100%	Feed with a minimum O ₂ flow of 40 L/min		

Performance table subject to changes. Always refer to product labelling.

Code	Description	Size	Box Qty.
CV204S	Ventumask 30 CPAP mask with integral Venturi flow driver and adjustable PEEP valve	Small	5
CV204M-L	Ventumask 30 CPAP mask with integral Venturi flow driver and adjustable PEEP valve	Medium / Large	5
CV204XL	Ventumask 30 CPAP mask with integral Venturi flow driver and adjustable PEEP valve	Extra Large	5
STS.15/30	Dual oxygen flow meter 15 + 30 L/min with O ₂ coupling		1



www.intersurgical.com/info/starmed

Video available at www.intersurgical.com

Make an enquiry

Dual flow meter manufactured by FLOW METER S.p.A.



Intersurgical Ltd, Crane House, Molly Millars Lane, Wokingham, Berkshire, RG41 2RZ, UK
 T: +44 (0)118 965 6300 F: +44 (0)118 965 6356 info@intersurgical.com www.intersurgical.com



Manufacturer:
Intersurgical S.p.A.
Mirandola (MO) Italy

Please think before you print
Save energy and paper. If you must print this information sheet please print it double sided.

IS1.43_INT • Issue 6 06.20